


Tube Plate 600 Evidence Sheet

—The effect of a plasticizing agent on metabolic stability assay—

Plates made from materials with inhibitors


Plates made from materials without inhibitor

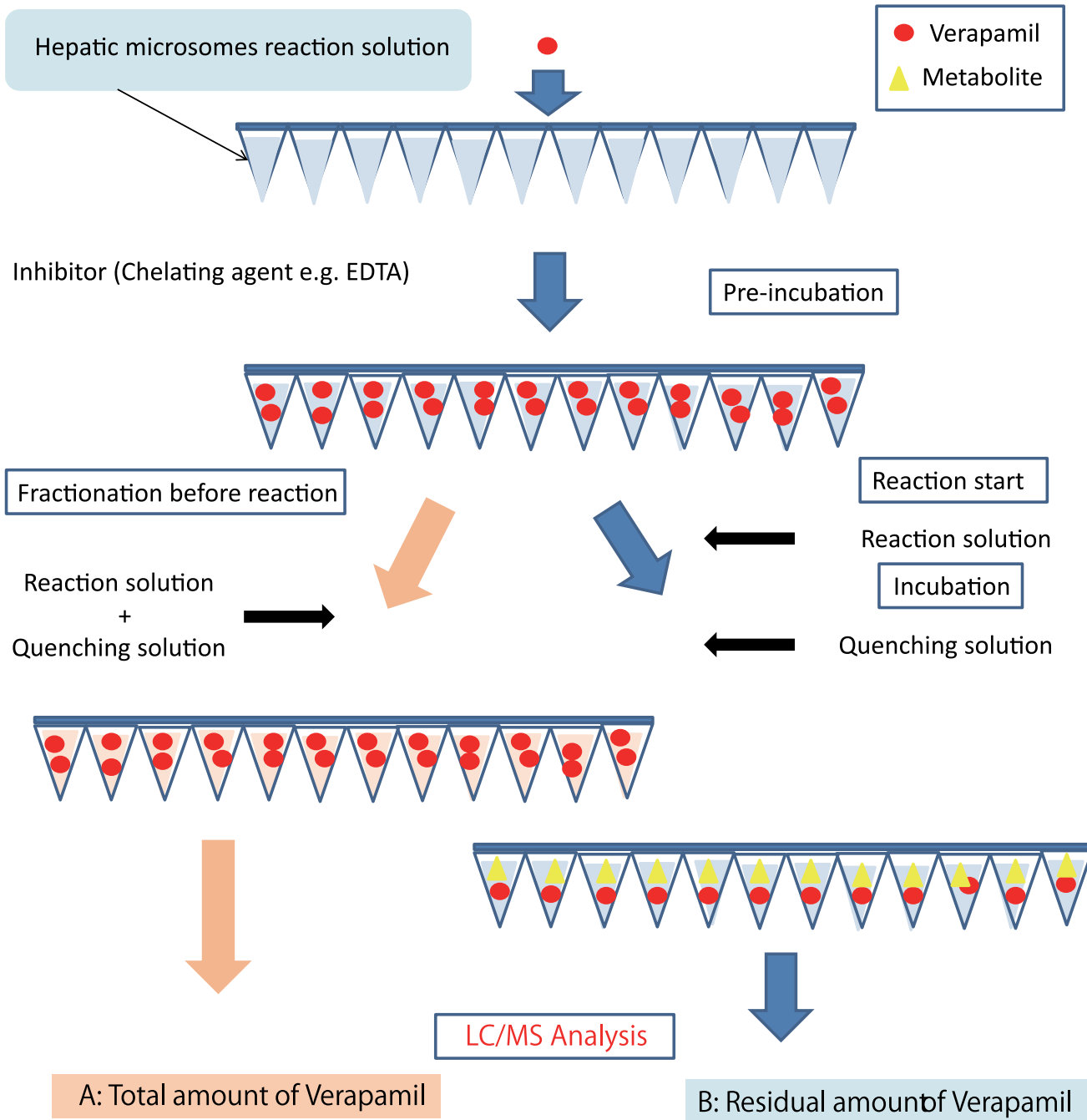

$$\text{Residual Rate} = \frac{B}{A} \times 100$$

Schematic drawing of an assay plate

	1	2	3	4	5	6	7	8	9	10	11	12
A												
B												
C												
D												
E												
F												
G												
H												

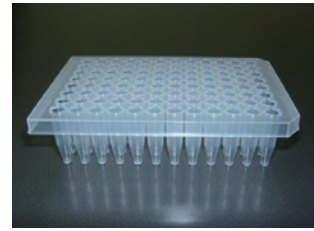

Plates made from materials with inhibitor

VS


Plates made from materials without inhibitor

Results!




【Result】

Plates made from materials with inhibitor

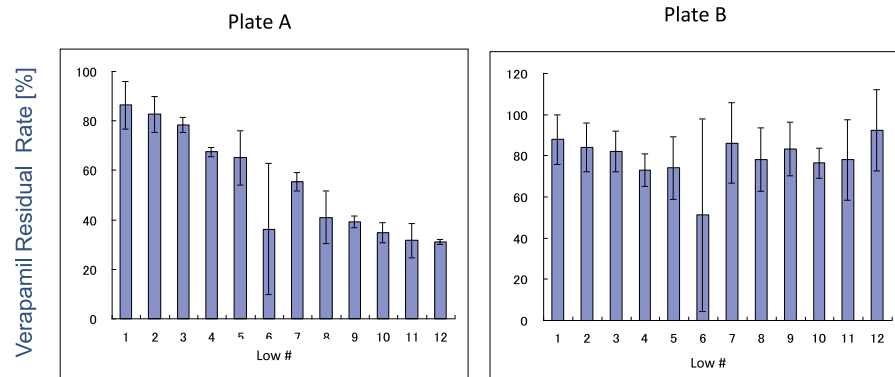

Data dispersion
between
wells & plates

Plates made from materials without inhibitor

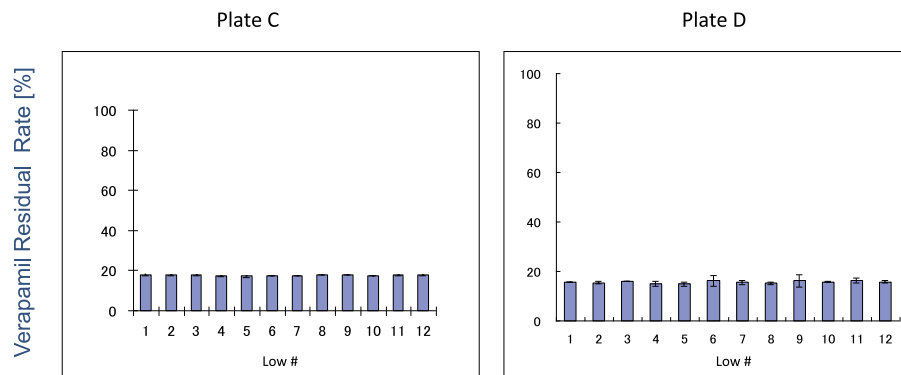

No data dispersion
between
wells & plates

Superior
Reproducibility